



CARE PUMP

## MOVE4 Sport line


<b>Number of chambers</b>	4
<b>Max. pressure [mmHg]</b>	220
<b>Number of proprietary therapeutic pre-programs</b>	6
<b>Smart gradient' function</b>	(mmHg, %, final pressure)
<b>Precise regulation of pressure and time</b>	in steps of 1, 5 or 10
<b>Setting of individual pressure in each chamber</b>	included
<b>Number of independent channels</b>	1
<b>Adjustment of the speed of the treatment</b>	1-5
<b>Cuffs</b>	overlapping system
<b>Display</b>	7" TN lcd color touch screen
<b>Weight [kg]</b>	2,5
<b>Dimensions [mm]</b>	290x130x130
<b>Power supply</b>	battery 21VDC/1A
<b>Estimated battery life from full charge</b>	4 h*
<b>Warranty</b>	1 year

\*test carried out using the "Warm-up, prevention" program and the pressure of 100 mmHg

# Accessories




## Included:

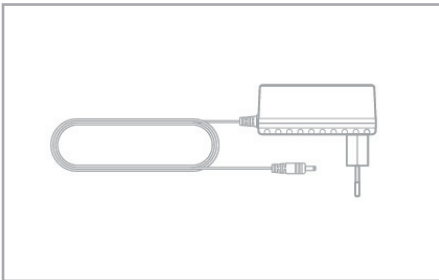


power supply

## Optional:

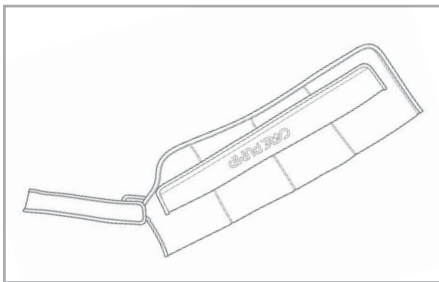


4-chamber arm cuff

CPARM4T

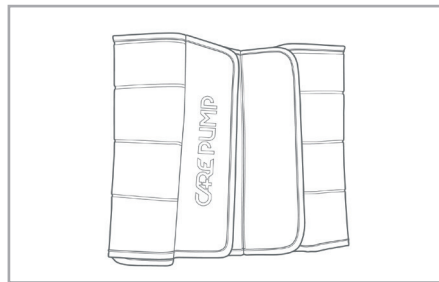


4-chamber waist cuff

CPWAIST4T

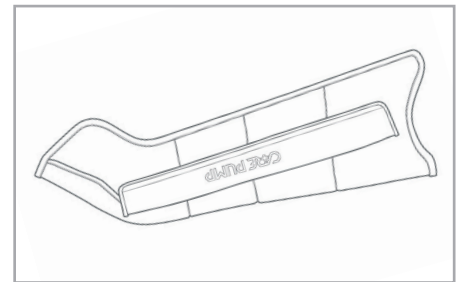


4-chamber left leg cuff

CPLEG4TL

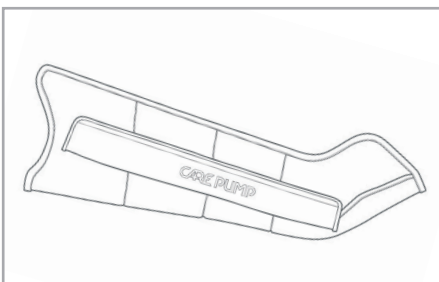


4-chamber right leg cuff

CPLEG4TR

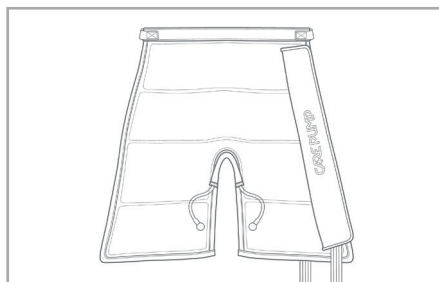


4-chamber short pants

CPSHORTPANTS4T

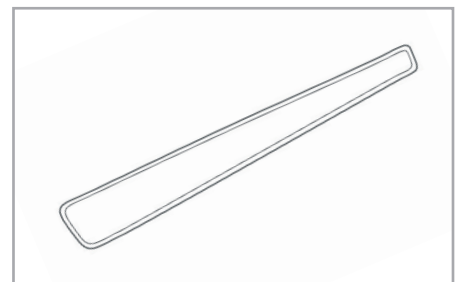


Extension zipper

CPLEGEXT4

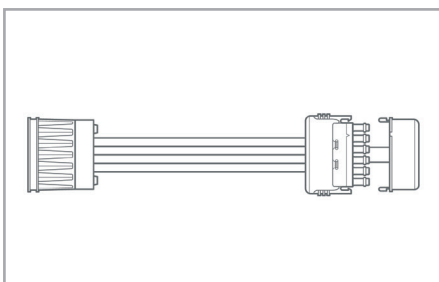


Double splitter

CPDOUBLECONNECTOR4

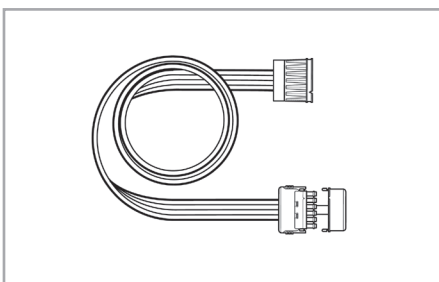

Double splitter (long - 150 cm)  
for 4-chamber devices

CPDOUBLECONNECTOR4L

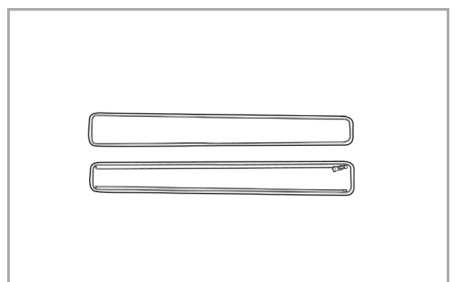

2 x Extension zippers  
for 4-chamber short pants

CPSHORTPANTSEXT4

# Therapeutic pre-programs



warm-up, prevention

- 
- 
- 
- 

Gradient	0 mmHg
Time	60 min.
Hold	2 s
Interval	2 s

regeneration, wellness

- 
- 
- 
- 

Gradient	4 mmHg
Time	45 min.
Hold	2 s
Interval	2 s

post-injury therapy

- 
- 
- 
- 

Gradient	4 mmHg
Time	30 min.
Hold	2 s
Interval	2 s

oedema therapy

- 
- 
- 
- 

Gradient	1 mmHg
Time	50 min.
Hold	3 s
Interval	3 s

all-chamber massage

- 
- 
- 
- 

Gradient	0 mmHg
Time	30 min.
Hold	2 s
Interval	3 s

mixed mode

**Possibility of combining any therapeutic pre-programs and selecting the number of repetitions for each of them.**